

Dental 3D Printers

Professional 3D printing technology
for dentists and dental technicians



CATALOG
2026

Asiga Max 2

Asiga Ultra

Form 4/4B

Form 4L/4BL

J3 DentaJet

J5 DentaJet

Dentacom Srl's 6 3D Printers

Printer	Technology	Pixel (µm)	Print volume	Materials	Strength
Asiga Max 2	DLP + SPS™	62 (shifting)	119×67×76 mm	>500 (open)	DLP precision + open system
Asiga Ultra	DLP 4K + SPS™	50 (shifting)	189×106×130 mm	>500 (open)	Large volume + 4K pixel shifting
Formlabs Form 4/4B	MSLA (LFD)	50	200×125×210 mm	10 + open	98.7% reliability — 100 mm/h
Formlabs Form 4L/4BL	MSLA (LFD)	46	353×196×350 mm	10 + open	Large format — high-volume labs
Stratasys J3	PolyJet	≤18 µm/layer	140×200×190 mm	3 simultaneous	Multi-material — final quality
Stratasys J5	PolyJet	≤18 µm/layer	140×200×190 mm	5 simultaneous	TrueDent® monolithic dentures

✓ All printers meet the requirements of the 2026–2028 Hyper-Depreciation scheme — Annex A, Law 232/2016

Key Asiga Technologies — What Really Sets Them Apart

1 SPS™

Smart Positioning System — Z-axis control

- Physical sensors measure the platform's actual position before every UV exposure
- Automatic correction of every deviation: no accumulated error on the Z axis
- Ultra: 4 proximity sensors for continuous self-calibration
- The Z height matches the CAD design exactly (tolerances $\pm 30 \mu\text{m}$)

2 Real resolution

Physical, shifted, or calculated pixel: not the same thing

- PHYSICAL pixel: every micromirror corresponds to a real pixel — geometric precision guaranteed. Used by Asiga.
- SHIFTED pixel (4K): the panel physically vibrates — doubling the real projected pixels. Used by Asiga.
- INTERPOLATED pixel: the software generates virtual pixels — estimated precision, NEVER used by Asiga
- A “30 μm calculated” can have the real sharpness of a “50 μm physical”

3 UV Radiometer

Measures actual power — not the nominal value

- The projector loses power over time: without compensation, the UV dose drops by 20-30% within a few months
- The radiometer measures the actual UV power after the entire optical chain (lens + glass + fluid)
- The system automatically adjusts exposure time to keep the dose constant
- Mechanical properties and part color remain stable for years

Asiga Max 2

Compact high-performance DLP — precision and open system

Print volume (X×Y×Z)	119 × 67 × 76 mm
Pixel size	62 µm (with 4K optical shifting)
Technology	DLP + SPS™ + UV radiometer
UV source	LED 385 nm
Crown & Bridge Kit	49×27×76 mm — tolerances ±30 µm
Materials	Open system — >500 certified resins
Connectivity	Wi-Fi, WirelessDirect, Ethernet
Dimensions / Weight	260×385×370 mm 19 kg

SPS™ — Z-axis control

Physical sensors verify the platform's actual position before every UV exposure. No estimation, no accumulated error.

Internal UV radiometer

Measures actual UV power at every layer. Automatically adjusts exposure: consistent quality for years. Quality never drops — only the print time lengthens slightly to compensate for the decrease in UV power.



Asiga Ultra

DLP 4K + SPS™ — large print volume, absolute precision


Print volume (X×Y×Z)	189 × 106 × 130 mm
Native physical pixel	25 µm
With 4K shifting	50 µm (physical pixels, not calculated)
Technology	DLP 4K + SPS™ + UV radiometer
UV source	LED 385 nm
Self-calibration	4 proximity sensors — automatic continuous calibration
Lid opening	Motorized — opens without touching it
Materials	>500 resins — open system

SPS™ + UV Radiometer

Z-axis control sensors + real UV measurement at every layer. Consistent quality for years.

Automatic calibration and opening

4 sensors ensure continuous calibration with no intervention. The lid opens automatically, hands-free.

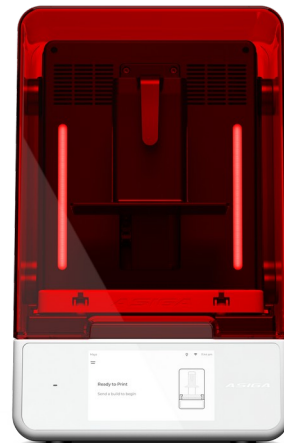

Asiga Crown Kit

Dedicated tray for crowns and bridges — prints in under 12 minutes

Complete kit for precision dental restorations

The Asiga Crown Kit is an interchangeable system (tray + build platform) designed to print crowns, bridges, inlays, onlays and veneers with maximum efficiency and minimal resin waste (only 30 ml per cycle). Compatible with Asiga Max 2 and Ultra.

⚡ Ultra-fast printing

Crowns and bridges completed in under 12 minutes per print cycle.

💧 Minimal resin consumption

Reduced tray volume of just 30 ml → lower waste and cost per unit produced.

🔧 Interchangeable

Mounts and removes in seconds on Max 2 and Ultra, no special tools required.

✅ Complete kit

Includes a dedicated build platform and a tray optimized for restorations.

🦷 Clinical applications

Crowns, bridges, inlays/onlays, veneers, esthetic mock-ups and temporary restorations.



Ultra Crown Kit Set



2L Ultra UltraGLOSS™
Build Tray



Max Crown Kit Set

Asiga UltraGLOSS™/Lift UltraGLOSS™ Build Tray

Pre-polished surfaces — zero manual post-processing

UltraGLOSS™ is a patented, ultra-reflective build tray. Glossy, pre-polished final result with no extra manual steps.

LIFT™ Technology

A liquid layer between the resin and the FEP film reduces peel forces. Less stress on parts, fewer surface defects.

Direct glossy finish

The surface comes out already polished. No need for abrasive steps or hand-polishing.

Faster printing

Reduced separation distances = faster cycles compared to standard trays.

Longer tray life

The liquid layer cushions mechanical stress: the tray lasts significantly longer.

Ideal for bite guards & splints

Recommended for bite guards, splints and clear resin orthodontic appliances.

Availability

Max / Max 2: 1 L • Ultra / Pro 4K: 2 L. Compatible with all Asiga resins.



Key Formlabs Technologies — What Really Sets Them Apart

1 LFS™

Low Force Stereolithography — near-zero peel force

- Patented flexible film: eliminates the peel force that causes distortion and print failures
- Result: isotropic surfaces, precise parts, 98.7% reliability across millions of prints
- Suited to delicate geometries: undercuts, thin bridges, full-arch crowns
- Less stress on the resin = more stable, repeatable mechanical properties

2 Integrated ecosystem

Hardware + software + materials developed together

- PreForm: free software with optimized automatic supports and print simulation
- Form Wash & Cure: dedicated post-processing with preset parameters for every resin
- Over 40 certified resins (biocompatible, models, temporaries, surgical guides)
- Dashboard fleet management to monitor multiple networked printers

3 Real speed

Up to 100 mm/h with no compromise on quality

- Form 4 and Form 4B print at 100 mm/h thanks to the high-speed LFS engine
- Full arch in ~17 minutes, single crown in ~8 minutes
- Z resolution: selectable 25–300 µm; XY: 50 µm (Form 4B)
- Print volume 20×18×24 cm: multiple parts in a single build

Formlabs Form 4/4B

Low-peel-force MSLA — certified speed and reliability

Print volume (X×Y×Z)	200 × 125 × 210 mm
Pixel size	50 µm
Technology	MSLA — Low Force Display (LFD)
Speed	Max 100 mm/h typical 40 mm/h
Materials	10 certified + open system
Reliability	98.7% — third-party certified
Connectivity	Wi-Fi, Ethernet, USB
Dimensions / Weight	418×372×620 mm 18,3 kg

Precision verified with GOM scanning

99.7% of surfaces within 100 µm of the original CAD, verified by comparing a scan of the printed model with the design file.

11 aligner models in 8 minutes

With dental aligner resins: 11 models in 8 minutes. Production capacity for a high-volume professional lab.



Formlabs Form 4L/4BL

Large format — 2.7× the volume of the Form 4/4B for high-volume work

Print volume (X×Y×Z)	353 × 196 × 350 mm (24.2 liters)
Pixel size	46 µm
Technology	MSLA — low peel force, 145 LEDs
Speed	Max 80 mm/h typical 24 mm/h
Materials	10 certified + open system
Precision	96.7% of surfaces within 50 µm of the design
Connectivity	Wi-Fi, Ethernet, USB
Dimensions / Weight	633×481×1020 mm 58,5 kg

2.7× greater volume

353×196×350 mm — 24.2 liters, the largest print volume in the dental MSLA category. Up to 75 models per cycle.



Same 4B platform

Same PreForm software, same materials, same workflow. No new training needed for anyone already using the Form 4/4B.

Key Strategies Technologies — What Really Sets Them Apart

1 PolyJet™

Simultaneous multi-material — up to 5 materials in a single print

- Jetting of liquid photopolymers with inkjet heads: layers as thin as $\leq 18 \mu\text{m}$, smooth as-printed surfaces
- 3–5 different materials for color, hardness and transparency in a single job
- TrueDent®: polychromatic with 500,000+ combinations for hyper-realistic esthetic dentures
- Soluble supports: complex geometries with no mechanical post-processing

2 Productivity

41 units/day — industrial-scale volume for the lab

- J5 DentaJet: up to 41 units/day, automatic mixing of different jobs on a single tray
- No manual material changes: multiple cartridges loaded, automatic switching
- GrabCAD Print: optimized nesting, production reporting, fleet management
- Measurable ROI: less labor, less waste, consistent output shift after shift

3 Biocompatibility

Materials certified for direct clinical use

- VeroDentPlus™, MED610™, TrueDent®: ISO 10993 certified and/or MDR Class IIa
- Stable colors, resistant to oral-cavity chemical agents
- Ultra-precise parts: controlled transparency and hardness
- Clinically validated: suited for dentures, surgical guides and long-term temporaries

Stratasys J3 DentaJet

Multi-material PolyJet — final quality straight off the printer

Print volume (X×Y×Z)	140 × 200 × 190 mm
Layer thickness	≤18 µm
Technology	PolyJet — UV jetting deposition
Simultaneous materials	3 materials at once
Main applications	Surgical guides, models, gingiva masks
Software	GrabCAD Print
Connectivity	Ethernet, Wi-Fi
Dimensions / Weight	620×620×590 mm 98 kg

3 materials at once

Surgical guides, models and gingiva masks are printed and cured in a single print cycle.

18 µm layers

Jetting deposition with 18 µm layers. Smooth surfaces, fine detail — no additional finishing needed.



Stratasys J5 DentaJet

5 materials — TrueDent® polychromatic monolithic dentures

Print volume (X×Y×Z)	140 × 200 × 190 mm
Layer thickness	≤18 µm
Technology	PolyJet — 5 materials at once
Cases per day	Up to 41 dental cases/day
Main materials	TrueDent®
Applications	Complete and partial removable dentures, temporaries
Software	GrabCAD Print
Dimensions / Weight	620×620×590 mm 228 kg

TrueDent® — Monolithic dentures

Polychromatic material that reproduces the natural shades of teeth. Monolithic dentures with no assembly required.



41 dentures per day

High-volume production for labs running multiple shifts. 5 materials simultaneously, no assembly required.

TrueDent — Like milled, better than printed

[Wear and durability comparable to milled restorations](#)

Wear and durability comparable to traditional milled teeth.

[Long-lasting color stability](#)

No premature yellowing: reliable esthetics even with coffee and wine.

[Fewer adjustments, greater satisfaction](#)

85% of TrueDent dentures need no adjustments, versus 14% of traditional dentures.

A maximum of 2 fitting appointments, versus 6 for traditional dentures.



Tax Benefits — 2026–2028 Hyper-Depreciation

180%

Uplift
maximum deductible

Gen. 2026

Effective
date

30 Set 2028

Booking deadline
(20% deposit)

L. 199/2025

Regulatory reference
Annex A, Law 232/2016

Key conditions

- Uplift on the tax-deductible cost (not an F24 tax credit)
- 3D printers listed in Annex A — Industry 4.0 technological capital goods
- Requirement: bidirectional connection to the company's IT system
- Met via Wi-Fi/Ethernet + Asiga Composer / PreForm / GrabCAD
- 5 mandatory filings to GSE required under the 2026 regulation
- Valid for companies and sole proprietorships carrying out an economic activity

Practical example (24% corporate tax)

Printer	Cost	Deduction ×180%	Corporate tax savings
Asiga Ultra (e.g.)	€ 16.500	€ 29.700	€ 3.168
Form 4/4B (e.g.)	€ 8.400	€ 15.120	€ 1.613

Figures are purely indicative. Actual savings depend on the tax rate (24% corporate tax for companies; personal income tax for sole proprietorships) and legal structure. Consult your accountant.

Why act by September 30, 2028: investments must be completed or booked (deposit ≥20%) by that date. The printers in this catalog already meet all the technical interconnection requirements.

Post-Processing Formlabs, Asiga & Stratasys

Asiga Cure • Form Cure 2nd Gen • Form Wash 2nd Gen • 5 L Pump-Fed Resin Dispensing System • Stratasys ProLight™

Asiga Cure

UV + Heat Post-Curing — Good Design Award 2024



Available sizes	1,5 L (194×139×65 mm) 2,5 L (194×139×132 mm)
UV source	LED 405 nm — high power, 360° immersion
Heating	Infrared — up to 80 °C, precise control
Curing atmosphere	Oxygen-free: vacuum or nitrogen
Connectivity	Wi-Fi, Wireless Direct, Ethernet/LAN
Programs	Automatic cycles for every resin — browser-based management
Dimensions / Weight	380×380×240 mm 15,6 kg
Design Award	Good Design Award Winner 2024

Oxygen-free Vacuum or nitrogen eliminate oxygen for maximum curing depth.	LEDs from all sides, 360° High-power LEDs distributed on every side ensure uniform curing on any geometry.	Automatic cycles per resin The system automatically selects the optimal parameters for every Asiga resin.	Fully open system Compatible with any photopolymer resin — no brand lock-in.
---	--	---	--

Formlabs Form Cure 2nd Gen

UV + Heat Post-Curing — integrated into the Formlabs ecosystem



UV source	LED 405 nm with rotating platform
Max temperature	80 °C (176 °F)
Compatibility	All Formlabs dental resins
Integration	Formlabs Dashboard — cloud monitoring
Curing protocols	Preset for every Formlabs resin
Recommended use	Form Wash + Form Cure — certified workflow
Availability	Included in the dental ecosystem package
Note	Resin-specific parameters in the IFU sheets

Certified cycles for every resin

Every Formlabs resin has preset times and temperatures, clinically validated for its MDR class.

Washing + curing

Form Wash cleans with isopropyl alcohol, Form Cure 2nd Gen completes curing — an integrated, traceable process.

Up to 80 °C

Heat accelerates conversion of residual material: mechanical properties at their peak.

Formlabs Form Wash 2nd Gen

Automatic washing — preset cycles for every dental resin



Type	Automatic wash station with agitation
Tank volume	4.5 liters
Solvent	Isopropyl alcohol (IPA) or Form Wash
Compatibility	Form 4/4B, Form 4L/4BL and the entire Form line
Cycles	Automatic — optimal duration for every resin
Basket	Dual: removable top + fixed bottom
Connectivity	Wi-Fi, Ethernet — integrated with the Dashboard
Recommended use	Form Wash → Form Cure 2nd Gen — full workflow

Automatic cycles per resin

PreForm and the Formlabs Dashboard automatically set wash duration for every resin — no risk of under- or over-washing.

Dual basket


The removable top basket and fixed bottom basket let you wash parts of different sizes in a single cycle.

MDR-certified workflow

Form Wash 2nd Gen + Form Cure 2nd Gen: post-production process validated for Class IIa dental resins.

Formlabs — 5 L Pump-Fed Resin Dispensing System

5 L jugs with pump — high volume, fewer interruptions

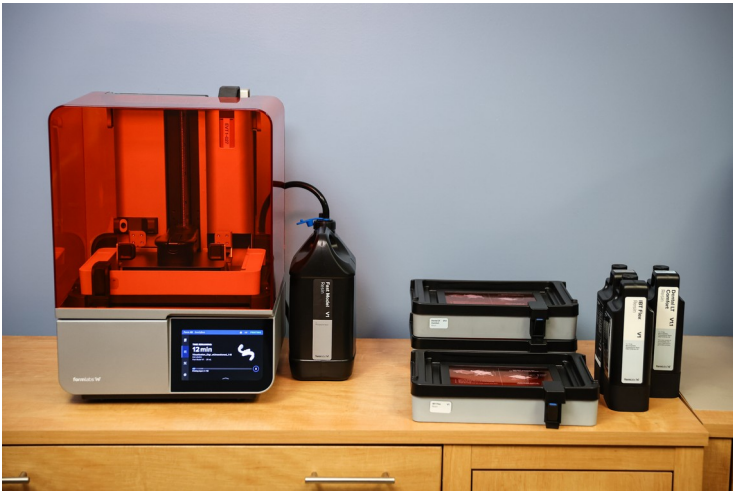

Container format	5-liter jugs (vs. 1 L cartridges)
System	Automatic dispensing pump
Compatibility	Form 4/4B, Form 4L/4BL — high-volume labs
Available resins	Fast Model Resin V1 and the dental resin line
Cost advantage	Lower cost per mL vs. single cartridges
Traceability	Resin level monitored by the Dashboard
Less waste	Less packaging and fewer production interruptions
Ideal for	Labs with high daily print volumes

Lower cost per mL

The 5 L jugs lower the unit cost of resin compared to 1 L cartridges — a concrete advantage for high-volume labs.

Continuous production

With 5 L of resin on hand, the printer keeps running without stopping to change cartridges. Fewer interruptions, more parts per day.

Simplified management

The Formlabs Dashboard monitors the level in the jugs and warns before they run out. Stock always under control.

Stratasys ProLight™ Curing Box

Post-curing for J3 and J5 DentaJet — large format



UV source	45 multi-directional 450 nm LEDs — 70 W electrical / 36 W radiant
Heating	2 × 500 W — fino a 80 °C (176 °F)
Curing volume	Cylinder ø 39.5 cm × h 32 cm
Rotating platform	Glass — uniform 360° exposure
Interface	LCD touchscreen
Connectivity	Wi-Fi 2,4 GHz, Ethernet 100 Mbit, USB 2.0
Dimensions / Weight	69 × 54 × 44,5 cm 24 kg
Operating temperature	18–28 °C

45 LEDs from every side

450 nm LEDs distributed throughout the chamber ensure uniform curing even for complex geometries.

Heat up to 80 °C

Two 500 W modules: precise temperature to complete PolyJet curing.


Rotating platform

The glass platform rotates during curing: no shadow spots, no partial curing.

Wi-Fi/LAN connection

Remote monitoring via GrabCAD Print — native integration with the Stratasys system.

Dental Resins

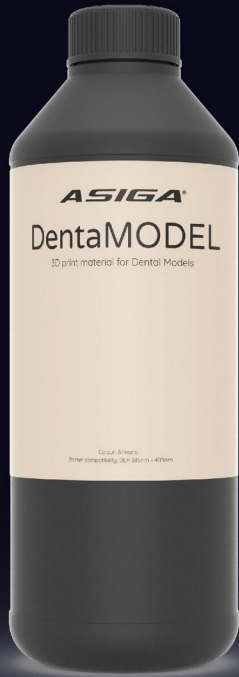


ASIGA

DentaMODEL

Study Model Resin

Ideal for printing anatomical study and working models with extremely high dimensional precision.

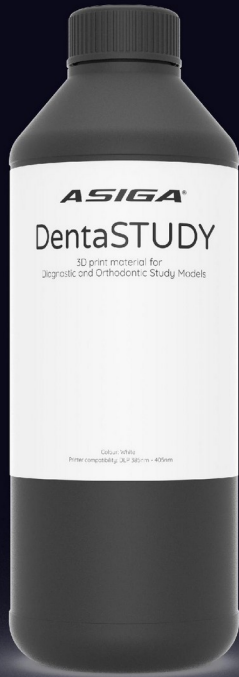

ASIGA

DentaSTUDY

Diagnostic Model Resin

Economical white resin for diagnostic and patient-communication models. Smooth surface, high detail and short print times for high productivity.




ASIGA

DentaFORM Beige

Models & Mock-ups (Beige)

Rigid resin for occlusal models, diagnostic mock-ups and esthetic prototypes. Beige shade that mimics the color of traditional plaster.

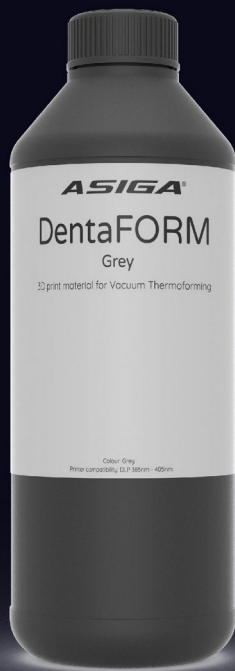


ASIGA

DentaFORM Grey

Models & Mock-ups (Grey)

Dark grey version of DentaFORM to improve visibility of occlusal detail during morphological analysis and implant planning.

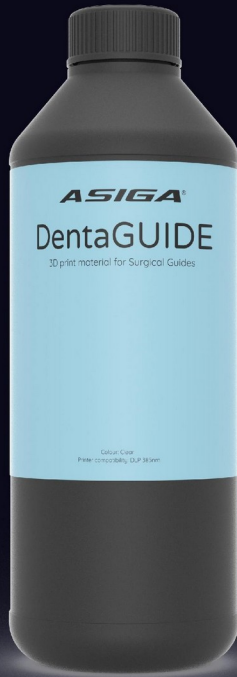


ASIGA

DentaGUIDE

Surgical Guides for Implants

Rigid, sterilizable resin for producing implant surgical guides. High mechanical strength during surgery, biocompatible for intra-oral use.

✓ Compliant with MDR 2017/745

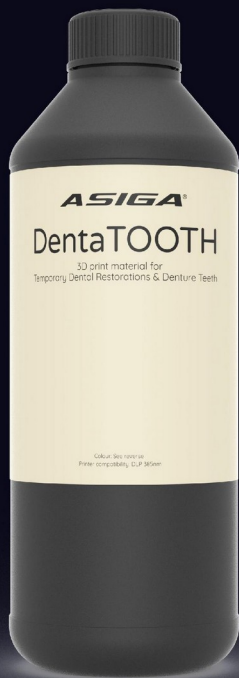


ASIGA

DentaTOOTH

Temporary Restorations

Resin for crowns, bridges and full-arch temporary restorations. Excellent flexural and abrasion resistance. Available in several VITA shades for immediate esthetic results.

✓ Compliant with MDR 2017/745

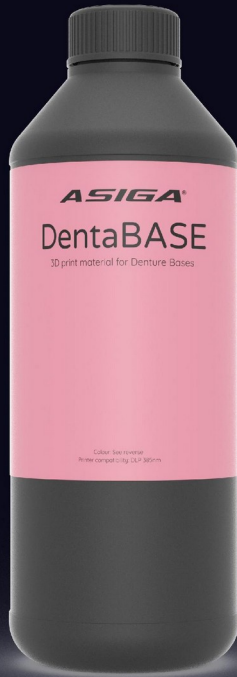


ASIGA

DentaBASE

Removable Denture Bases

High-strength rigid resin for printing complete and partial denture bases. Excellent surface finish, compatible with both resin and ceramic teeth.

✓ Compliant with MDR
2017/745

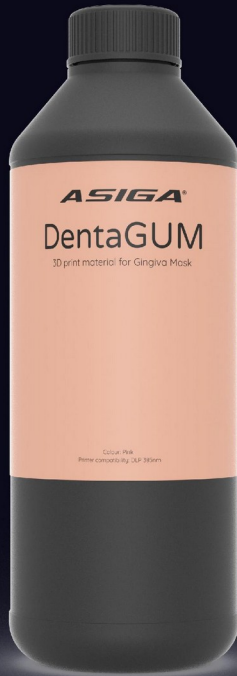


ASIGA

DentaGUM

Gingiva Simulation

Flexible, gingiva-colored resin for printing gingiva masks on hybrid models. Enables realistic soft-tissue simulation for implant planning.

✓ Compliant with MDR
2017/745

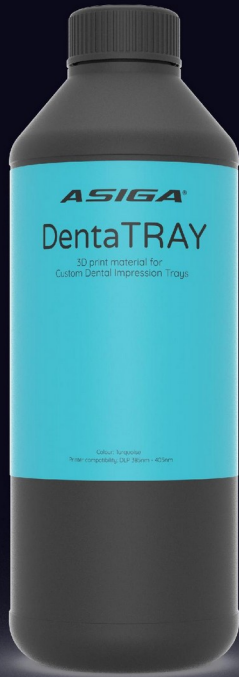


ASIGA

DentaTRAY

Custom Impression Trays

Clear resin for printing individual trays for precision impressions. Light, rigid and compatible with all major impression materials.

✓ Compliant with MDR
2017/745

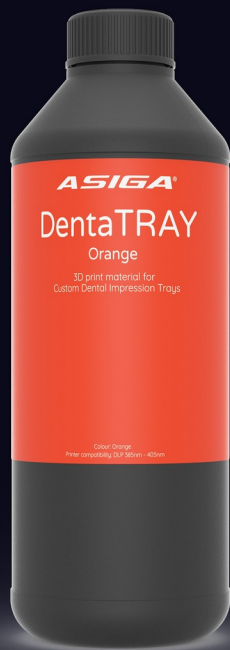


ASIGA

DentaTRAY Orange

Custom Impression Trays (Orange)

Fluorescent orange variant of DentaTRAY for greater visibility chairside and instant identification in the lab. Same mechanical properties as the standard version.

✓ Compliant with MDR 2017/745

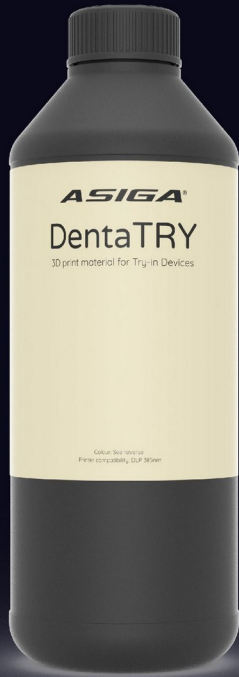


ASIGA

DentaTRY

Esthetic Try-in

Resin for printing provisional try-ins. Enables esthetic and functional testing before the final restoration, reducing chairtime.

✓ Compliant with MDR
2017/745

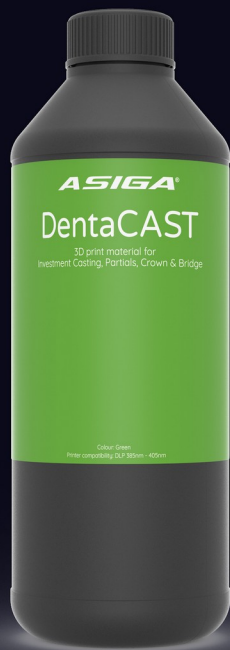


ASIGA

DentaCAST

Castable Restorations

Clean-burning castable resin for lost-wax casting. High dimensional precision with zero ash residue, ensuring clean castings.

✓ Compliant with MDR
2017/745

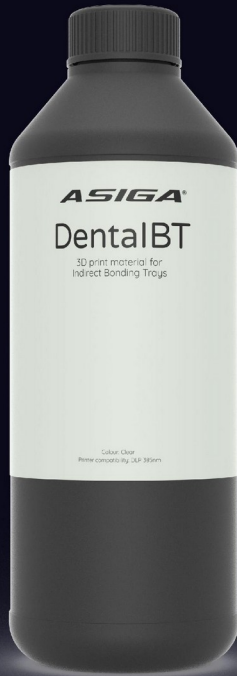

ASIGA

DentalBT

Indirect Bracket Bonding Trays

Rigid resin for producing indirect bonding trays in orthodontics. Precise bracket placement, fully digital workflow.

✓ Compliant with MDR 2017/745

FORMLABS

Precision Model Resin

High-precision working models

Resin for study and working models with high surface detail.

Compatible with: Formlabs Form 4/L



FORMLABS

Fast Model Resin

Diagnostic, orthodontic and aligner thermoforming models

Fast resin for diagnostic and orthodontic models. The fastest choice for high-volume labs.

Compatible with: Formlabs Form 4/L



FORMLABS

Flexible 80A Resin

Gingiva masks

Flexible resin for gingiva masks on hybrid models. Requires Magenta + White + Yellow pigments to achieve the gingival pink.

Compatible with: Formlabs Form 3/4/L



FORMLABS

Castable Wax Resin

Wax models for casting and pressing



Castable resin for lost-wax casting.
Vaporizes completely.

Compatible with: Formlabs Form 3/4/L

FORMLABS

Surgical Guide Resin

Implant surgical guides

Rigid, sterilizable resin for implant surgical guides. High mechanical strength during surgery. Class Ia certified.

✓ **Compliant with MDR 2017/745**

Compatible with: Formlabs Form 3B/4B/BL



FORMLABS

IBT Flex Resin



Trays for indirect bracket bonding and direct composite injection

Resin for producing flexible trays for indirect orthodontic bonding, also used for injecting composite restorations. Class Ia certified.

✓ **Compliant with MDR 2017/745**

Compatible with: Formlabs Form 3B/4B/BL

FORMLABS

Dental LT Clear v2



Rigid bite guards and occlusal splints

Rigid clear resin for bite guards and occlusal splints.
Class IIa certified.

✓ **Compliant with MDR 2017/745**

Compatible with: Formlabs Form 3B/4B/BL

FORMLABS

Dental LT Comfort



Flexible bite guards and occlusal splints

Flexible clear resin for bite guards and whitening trays.
Class IIa certified.

✓ **Compliant with MDR 2017/745**

Compatible with: Formlabs Form 3B/4B/BL

FORMLABS

Premium Teeth HT

Teeth for dentures and temporary restorations

Resin for teeth in complete dentures, temporary crowns/bridges and all-on-x immediate loading. Available in 4 shades (A2/A3/B1/Bleach).

✓ **Compliant with MDR 2017/745**

Compatible with: Formlabs Form 3B/4B/BL



DETAX

Freeprint Crown

Definitive crowns and bridges

Resin for definitive crowns and bridges with high translucency and color stability. Class IIa certified.

✓ **Compliant with MDR 2017/745**

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



DETAX

dx direct aligner

Directly printed aligners

Resin for direct printing of aligners. Class IIa certified.

✓ **Compliant with MDR 2017/745**

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



DETAX

dx denture flex

Flexible partial dentures (definitive use)

Resin for long-term definitive flexible partial dentures.
Class IIa certified.

✓ **Compliant with MDR 2017/745**

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



VOCO

V-Print SG

Implant surgical guides

Resin for implant surgical guides. Sterilizable and biocompatible, Class Ia certified.

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



VOCO

Conforme
MDR 2017/745

V-Print splint 2.0

Rigid bite guards and occlusal splints

Rigid resin for rigid bite guards and occlusal splints.
Class IIa certified.

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



VOCO

Conforme
MDR 2017/745

V-Print splint comfort

Flexible bite guards and occlusal splints

Flexible resin for bite guards with superior comfort.
Elastic properties for greater joint tolerance. Class IIa
certified resin.

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



VOCO

Conforme
MDR 2017/745

V-Print c&b temp

Temporary crowns and bridges

Resin for temporary crowns and bridges. Complements options requiring additional shades or higher mechanical specifications. Class IIa certified resin.

Compatible with: Asiga MAX UV and Ultra • Formlabs Form 4B



VOCO

Conforme
MDR 2017/745

V-Print Try-In

Resin for esthetic try-ins

Rigid resin for esthetic try-ins and test specimens for partial or complete dentures. Class Ia certified.

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



VOCO

Conforme
MDR 2017/745

V-Print dentbase

Complete denture bases

Rigid resin for complete denture bases. Class IIa certified.

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



BEGO

Conforme
MDR 2017/745

VarseoSmile Temp

Temporary crowns and bridges

Resin for temporary crowns and bridges. Complements labs that need additional shades or higher mechanical specifications than DentaTOOTH and Premium Teeth HT. Class IIa certified resin.

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



BEGO

VarseoSmile Crown plus

Definitive crowns and bridges

Resin for definitive crowns and bridges. Class IIa certified.

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



BEGO

Conforme
MDR 2017/745

VarseoSmile TriniQ

Definitive bridges up to 3 units

Resin for definitive crowns and bridges up to 3 units.
Class IIa certified.

Compatible with: Asiga MAX and Ultra • Formlabs Form 4B



NEXT STEPS

Let's talk about your digital project

EMAIL: info@dentacom.it

WHAT WE OFFER

Personalized ROI analysis

Tailored quotation

Free sample print

Personalized pre-sale training

Dedicated after-sales support

Hands-on technical courses

Asiga Max 2
& Ultra

Form 4/4B
& 4BL

Stratasys J3
& J5

DENTACOM SRL

May 2026 catalog edition • Specifications subject to change without notice